
- A** Reservoir
- B** Max Fill Line
- C** Tube
- D** Button
- E** Pulp Housing
- F** Nozzle

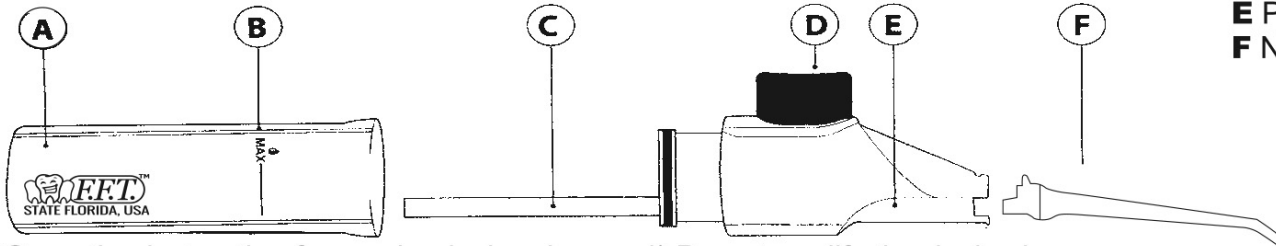Operating instruction for mechanical oral irrigator Individual Favorite Irrigator™ model F.F.T.-IFI-1000(F.F.T. Individual Favorite Irrigator 1000). Before using this product, please read this manual carefully and keep it for future reference.

Content

- 1) safety instructions
- 2) purpose and contraindications
- 3) completeness and technical characteristics
- 4) preparation of the device for operation
- 5) the use of the device
- 6) after using the device
- 7) rules of care, cleaning and disinfection
- 8) warranty obligations

Thank you for buying F.F.T. Individual Favorite Irrigator 1000 made by company F.F.T. (Favorite For Teeth), USA. This device provides effective cleaning of the space between the teeth and removal of plaque. Its use prevents the formation of Tartar and periodontal disease. F.F.T. Individual Favorite Irrigator 1000 is a necessary device for the care of dentures, crowns and braces. It can also be used to massage the gums. Daily home gum massage is an additional measure of caries prevention and the formation of inflammatory processes in the oral cavity.

THE SAFETY INSTRUCTIONS

ATTENTION! Do not throw the F.F.T. Individual Favorite Irrigator 1000 into a fire, do not heat, and do not leave the device in places with high air temperature. These actions can lead to a damage of the device, as well as reduce its life.

WARNINGS:

- 1) Do not use device if the tank, pump body, tube or nozzle of the irrigator is damaged.
- 2) Do not break, bend, or pull the irrigator nozzle with a force. This may reduce its life.
- 3) Keep the device out of the reach of children. It may result in accidents and dangerous situations, such as accidental swallowing of accessories or removable parts of the device by the child.

- 4) Do not modify the device in any way. Also do not disassemble or repair it yourself.
- 5) Do not store the device where it may fall and break.
- 6) Do not use the device if the nozzle is lost.

CAREFULLY:

- 1) Children under 14 years of age and persons with disabilities who cannot operate the device on their own, as well as persons with impaired oral sensitivity should not use this device. People undergoing oral treatment or experiencing pain in their teeth or gums may not always use an irrigator. Before using the device, you must contact your dentist.
- 2) Use the device only for cleaning the oral cavity. Water pressure directed at the eyes, nose, ears or throat can cause injury.
- 3) Do not press the irrigator nozzle against your teeth or gums during use, as this may cause damage.
- 4) Do not fill the tank with water whose temperature exceeds 40 °C.
- 5) Keep the device clean and use only fresh water or mouthwash of the same consistency as water to clean the mouth.
- 6) In case of non-compliance with sanitary and hygienic standards bacteria can develop in the device, which can be harmful to health. It can also cause an unpleasant odor.

APPOINTMENT AND CONTRAINDICATIONS

F.F.T. Individual Favorite Irrigator 1000 is produced for the prevention of bleeding gums, caries, periodontal diseases, gingivitis, medical irrigation and massage of the oral mucosa; effective cleaning of orthodontic systems, as well as for the care of the oral cavity and interdental spaces. It is produced for individual use only.

CONTRAINDICATIONS:

During the first few applications of the irrigator, there may be slight bleeding of the gums. If you have healthy gums, any bleeding of the gums should stop within 1-2 weeks. If bleeding continues after 2 weeks of use, discontinue use of the irrigator and consult your dentist.

Children age less than 14 years

Persons with disabilities

Persons suffering from neurological diseases

Persons suffering from mental illness

PACKAGE CONTENTS AND CHARACTERISTICS:

SPECIFICATIONS CONTENT:

- 1) the body of the irrigator
- 2) water tank
- 3) 1 standard nozzle
- 4) instruction manual
- 5) box

- 1) dimensions (W × d × h): 47 × 32 × 213 mm
- 2) weight: 94 g
- 3) water tank capacity: 50ml
- 4) pressure: from 45 to 75 psi


PREPARING THE DEVICE FOR USE
Before using the F.F.T. Individual Favorite Irrigator 1000, disconnect the water tank from the body and submerge the pump body, tube and nozzle completely in water. Press the irrigator button and hold them under water for at least 60 seconds until it starts to bubble. Then lift the end of the nozzle above the water and press the button to ensure the device is working properly and has a steady jet stream.

USE OF THE DEVICE

- 1) Holding the F.F.T. Individual Favorite Irrigator 1000 upright, gently turn the tank counter clockwise and separate it from the pump housing.
- 2) Fill the reservoir of F.F.T. Individual Favorite Irrigator 1000 up to the maximum filling line with either water or your choice mouthwash. Always use only fresh water or mouthwash with the same consistency as water.
- 3) Replace the pump housing by inserting it into the tank and gently turn clockwise to lock into place.
- 4) Using F.F.T. Individual Favorite Irrigator 1000 keep the device in a strictly vertical position, do not overturn the irrigator horizontally. Otherwise, the water will leak into the sections under the tank.
- 5) Place your thumb, as shown on photo (see above), on the button. Insert the nozzle of the irrigator into the mouth. To prevent splashing, cover your mouth slightly during irrigation. Press the button several times until the liquid comes out in a smooth flow. You can adjust the flow rate by changing the force you apply when you press the button. For best results, use a different pressure on the button and change the angle of the water jet by moving the F.F.T. Individual Favorite Irrigator 1000 over your mouth.
- 6) To clean the tooth gaps: place the nozzle in front of the teeth on the gum line, and press the button. For optimal results, use the appliance on both sides of the tooth, continuing to use the F.F.T. Individual Favorite Irrigator 1000 until the oral cavity is completely cleaned.
- 7) To clean bridges, braces, implants or other devices, place the nozzle under the bridge or other device as perpendicularly as possible and press the button. Move slowly along the tooth row. Special attention should be

! When operating F.F.T.-IFI-1000™

DO NOT tilt beyond horizontal otherwise water will leak from the slits under the reservoir

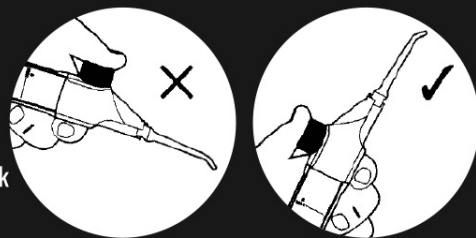paid to the treatment of interdental space and the area around the braces. For best results, start with the area of the molar(back teeth), moving gradually to the front teeth. Follow along the gums, stopping at the interdental spaces. Continue until all areas around and between the teeth are disinfected. However, remember that irrigation of the mouth will not replace brushing your teeth.

AFTER USING THE DEVICE

- 1) Disconnect the tank and pour the water away. This will help prevent contamination and growth of bacteria in the device.
- 2) Press the button several times to pour away the remaining liquid from the nozzle.
- 3) Wipe the irrigator with a dry cloth.
- 4) Storage and operating conditions of the device:
-0°C- + 40 °C

CARE RULES, CLEANING AND DISINFECTION

When rinsing the F.F.T. Individual Favorite Irrigator 1000, use water with a temperature not exceeding 40°C. Never use caustic soda or any abrasive substances to clean the device, as they can damage it. If you do not intend to use the device for a week or longer, dry the container inside and out.

WARRANTY OBLIGATION

The manufacturer guarantees the warranty for a period of 12 months from the date of purchase of the device. If a manufacturing defect is found during the warranty period, the defective instrument will be repaired and, if repair is not possible, replaced free of charge. Guarantee obligations are valid only if there is a cash receipt.

Warranty service is not provided:

- 1) in the presence of traces of mechanical action on the body of the device, dents, cracks, chips, etc,
- 2) in the presence of traces of opening the body, traces of attempts to repair,
- 3) in the presence of traces of ingress into the body or exposure to aggressive media,
- 4) in other cases of violation by the consumer of the rules of storage, transportation and technical operation of the device provided by the rules stated in the instruction manual.


The warranty does not apply to the packaging of the device. The company reserves the right to make full or partial changes to the products without prior notice and in accordance with the production requirements.

Manufacturer:

Favorite For Teeth (F. F. T.), US-United States, 2283 NE 164th st North Miami Beach, Florida State, 33160.

Developer:

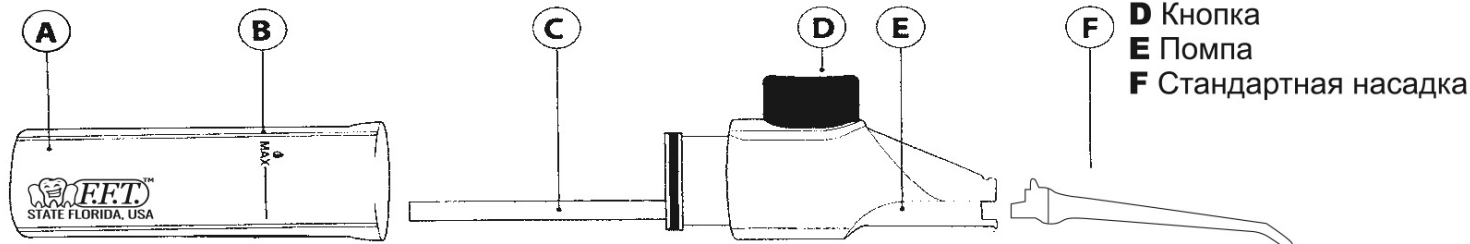Cron-USA LLC, US-United States, 2283 NE 164th st North Miami Beach, Florida State, 33160.

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